

AmbirScan for Athena User Guide



AmbirScan®



AmbirScan for athenahealth User Guide

Contents

- Software Installation Requirements 3
- Software Installation..... 3
- System Setup 4
- Opening and Running the Application..... 5
- QuickScans 6
- Selecting a QuickScan from Kiosk Mode..... 6
- Selecting a QuickScan from Interface 6
- Creating New QuickScans 7
 - Keeping an Existing QuickScan Setting 7
 - Creating a New QuickScan 7
- Deleting a QuickScan 8
- Scanning Workflows..... 9
 - Driver’s Licenses and Insurance Cards 9
 - Document Scanning 9
- Settings Menu 10
- QuickSettings Menus 13
 - Athena Card – Scanning 13
 - Athena Card – Output 14
 - Card – Imaging 15
 - Document – Output 17
- Troubleshooting and Technical Support 18



AmbirScan for athenahealth User Guide

Software Installation Requirements




You must have local administrative rights for initial installation of your software. If you are unsure of your security permissions, please contact your local IT professional.

Software Installation

To install your software, follow these steps:

1. **Download** the software from the link provided in your email
2. **Double click** AmbirScan for Athena Setup
3. Follow the prompts to complete installation
4. **Enter the License Key** provided in your email
5. **Enter the Unlock Code sent via SMS to the mobile phone numbers specified by your practice in the Business Associate Agreement with Ambir Technology. For security, Unlock Codes may only be sent via SMS to the specified phone numbers.**

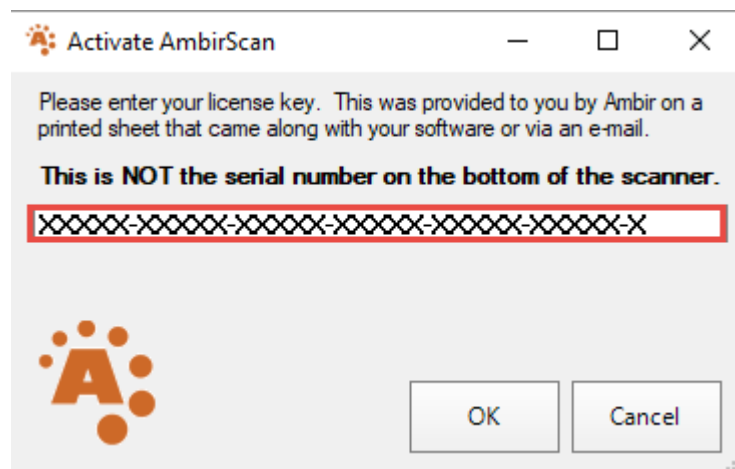

You may receive an error due to firewall settings. For assistance with error messages during installation, contact Ambir technical support



AmbirScan for athenahealth User Guide

System Setup

The below default folders are created by the AmbirScan for Athena software installer. The QuickScan settings in the AmbirScan for Athena application by default are setup to save cards and documents into folders on your local computer. The Card folder is where card images are stored, and the Athena Scans folder is for document scanning.

Affirm the following local folders are available:

- C:\Card
- C:\Athena Scans



If you do not see these folders, you must create them manually on your computer's C drive.



AmbirScan for athenahealth User Guide

Opening and Running the Application

The default setting for AmbirScan for Athena is to run in **Kiosk** mode. This means the software runs as a background process in your system tray. The AmbirScan for Athena icon will display in the system tray.

Note: Users will not see or interact with the interface unless they manually open the software.



AmbirScan for athenahealth User Guide

QuickScans


A **QuickScan** is a pre-set scanning profile (specifying file type, resolution, color, etc.) within the AmbirScan for Athena software. Users can edit, create, or delete QuickScans within the software.

The default selected QuickScan for AmbirScan for Athena is the Athena Card QuickScan. Users must change to a Document QuickScan when scanning any media other than cards.

Note: *The last QuickScan that was used will remain selected when AmbirScan for athena relaunches.*


Selecting a QuickScan from Kiosk Mode

1. **Right-click** on the AmbirScan for Athena icon in the system tray.
2. Select **QuickScans**, then choose the appropriate QuickScan for your scan. Once selected, a checkmark will appear to the left of the enabled QuickScan.

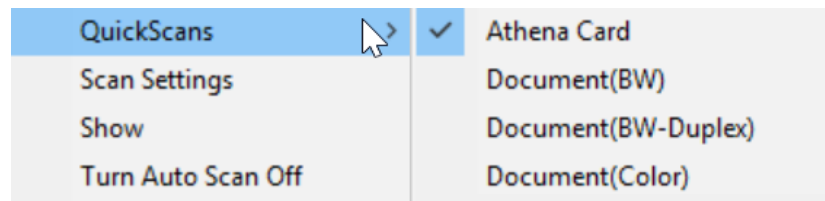

Selecting a QuickScan from Interface

1. **Right-click** on the AmbirScan for Athena icon in the system tray.
2. Select **Show**.
3. Select **QuickScans**, then choose the appropriate QuickScan for your scan. Once selected, a checkmark will appear to the left of the enabled QuickScan.
4. Click **Kiosk** to minimize the application again.

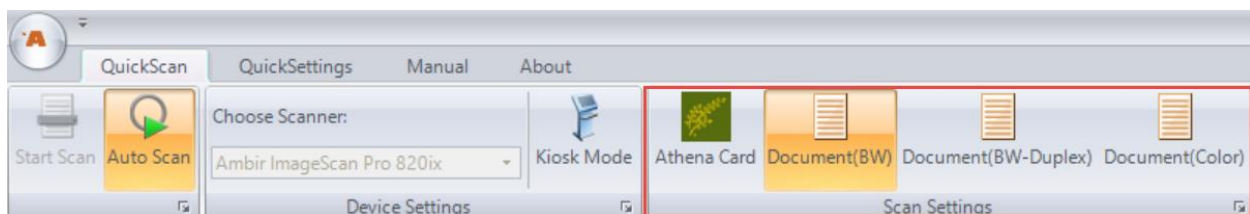

⚠ Clicking the “X” at the top of right of the window will exit the application completely. You will need to relaunch the software from your desktop to begin scanning again.



AmbirScan for athenahealth User Guide

Creating New QuickScans


There are two ways you can create a new QuickScan: editing an existing QuickScan and saving as new, or using the Wizard to start from the beginning.

Keeping an Existing QuickScan Setting

Users can create new QuickScans by editing an existing QuickScan and using the **Save as New QuickScan** feature in the **QuickSettings** menu.

To create a new QuickScan with AmbirScan for Athena settings, follow the below steps:


1. **Open** AmbirScan for Athena
2. **Click** the **QuickSettings** tab
3. **Select** an existing QuickScan
4. **Edit** the QuickScan with your changes
5. **Click Save as New QuickScan** in the **Operations** section
6. Your new QuickScan will appear in the **QuickScan** menu

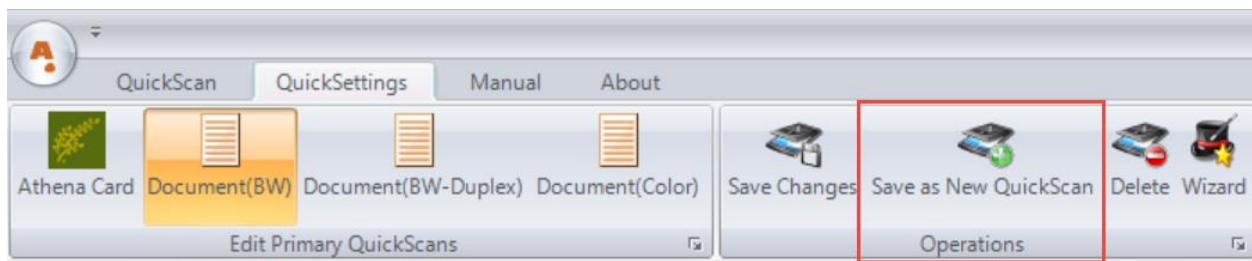

Creating a New QuickScan

To create a new QuickScan, follow the below steps:

1. **Open** AmbirScan for Athena
2. **Click** the **QuickSettings** tab
3. Click the **Wizard** function in the **Operations** section
4. **Name** and select and **Icon** for your QuickScan
5. **Select** your save path to the appropriate Athena folder on your C drive (*default save location is C:\Temp*)
6. Follow the prompts through the Wizard to select your settings
7. Click **Finish**
8. Your new QuickScan will appear in the **QuickScan** menu

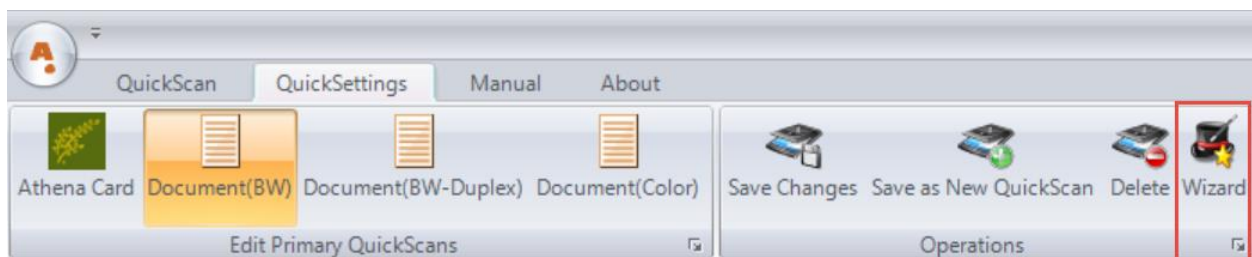




AmbirScan for athenahealth User Guide

Deleting a QuickScan

If you need to delete a QuickScan, you can do so by selecting the **Delete** button in the **QuickSettings** menu.

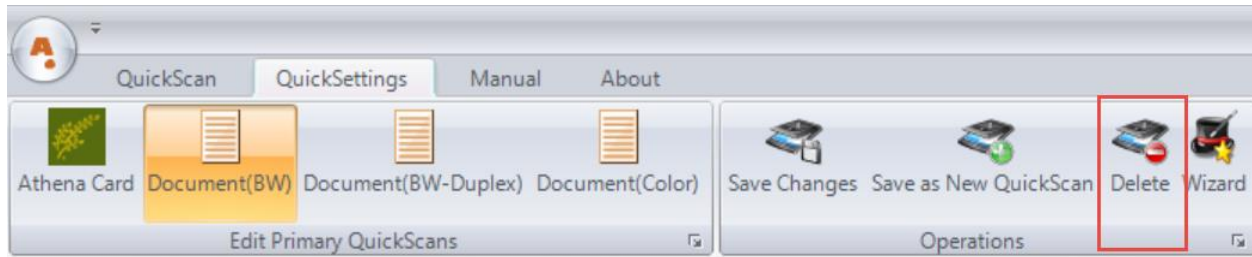

! Any changes made to QuickScans (Scanning, Output, Imaging,) will not take effect until you click on the Save Changes icon within the Operations box in the title bar area.



AmbirScan for athenahealth User Guide

Scanning Workflows

Driver's Licenses and Insurance Cards

To scan these supported cards, complete the following steps:

1. Confirm the **Athena Card QuickScan** is selected (*Verify by right-clicking the AmbirScan for Athena icon in the system tray and select QuickScans. A checkmark will be present next to the active QuickScan*).
2. **Insert** the ID card in the scanner. The software will automatically pull the card through the scanner
3. **Select** the correct patient from the search results
4. Click **Add Card Image** in patient chart under the correct section

Document Scanning

These are the default QuickScans for scanning documents via AmbirScan for Athena:

- Document(BW)
- Document(BW–Duplex)
- Document(Color)



AmbirScan for athenahealth User Guide


Settings Menu

The Settings menu is where specific administrative settings for the software are stored. The information below provides details as to what the settings are and how they are used.

There are two ways to access the Settings Menu:


From Kiosk Mode

1. **Right-click** on the AmbirScan for Athena icon in the system tray
2. Select **Scan Settings**

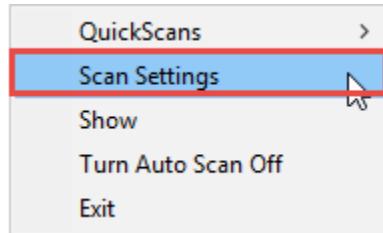

From AmbirScan for Athena Interface

1. **Right-click** on the AmbirScan for Athena icon in the system tray
2. Select **Show**
3. Click on the **AmbirScan icon** in the top left corner of the screen
4. Click **Settings**

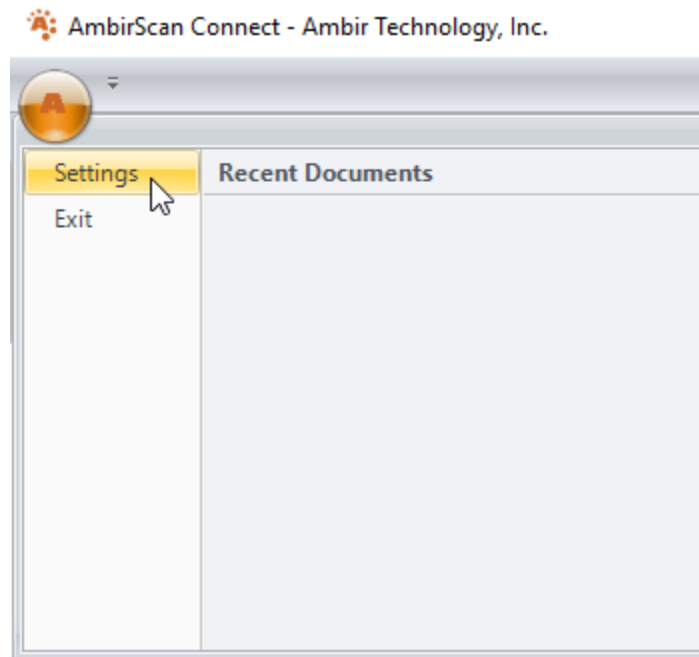



AmbirScan for athenahealth User Guide

Settings

The screenshot below displays the default settings in the **Settings** Menu:

Post Capture section

Show Document: When enabled, the document will open on your screen after scanning, in the application specified by Windows for the file type of the scanned document. This function is only available when using a QuickScan that it setup to save a local copy.

Auto Scan Section

Auto Scan (default): Any card/document inserted into the selected scanner will be pulled through automatically

Note: Do not adjust these settings unless directed by Ambir Support



AmbirScan for athenahealth User Guide

Administration Section

Start in Kiosk Mode (default): AmbirScan for Athena runs in the system tray

Auto Restart Kiosk Mode After User Login (default): Application will run in Kiosk mode after login

QuickScans Available to All Users (default): All users will have access to all QuickScans

QuickScans Available Only to Me: QuickScans created for use under specific profiles

Perform Logging: Only use this function when directed by Ambir Support

Debugging Console: Advanced functionality. Only for use when directed by Ambir Support

Calibrate: Calibration should be performed when image quality degrades. Click [here](#) for details on calibration.

Clean: Along with calibration, cleaning should be performed when image quality degrades. Click [here](#) for details on cleaning.

FTP Section

FTP Upload: A copy of the document will also be uploaded to a set FTP location

Printing Section

Print 1:1 (default): Print the actual size of the document

Print Fit (default): Scales to the size of the sheet

Show Print Dialog (default): Displays the print screen

Show File Name (default): Displays the file name

Email Section (this function is available when using a QuickScan that is setup to save a local copy)

SMTP Username: Enter your SMTP username

SMTP Password: Enter your SMTP password

SMTP Server Name: Enter the server name for your

Note: Only SMTP email configurations are supported.




AmbirScan for athenahealth User Guide

QuickSettings Menus

The **QuickSettings Tab** contains several sub-menus for settings for each QuickScan: **Scanning**, **Output**, **Imaging**. The following sections describe each of these menus in further detail for Card and Document QuickScans.

Athena Card – Scanning

The screenshot below displays the default settings for the **Card QuickScan Scanning Menu**:

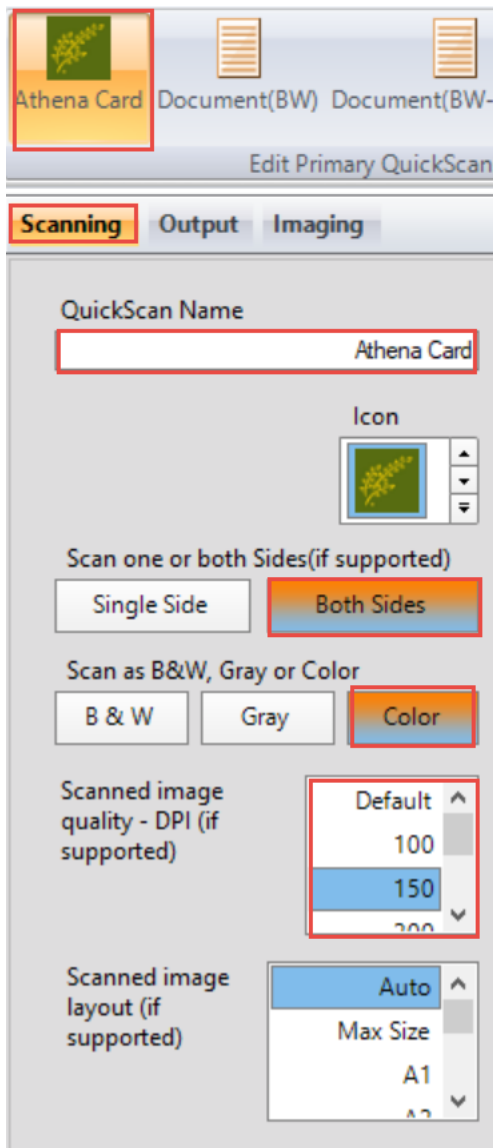

QuickScan Name: Name of QuickScan (Athena Card)

Both Sides: Selected for double sided scanning (default for Athena Card scans)

Scanned image quality – DPI: Image scan quality. The default for Card scans is 150 DPI.

Scanned image layout: The sizes of the media you are scanning. Auto is the default selection for Card scans.



AmbirScan for athenahealth User Guide

Athena Card – Output

The screenshot below displays the default settings for the **Athena Card QuickScan Output Menu**:

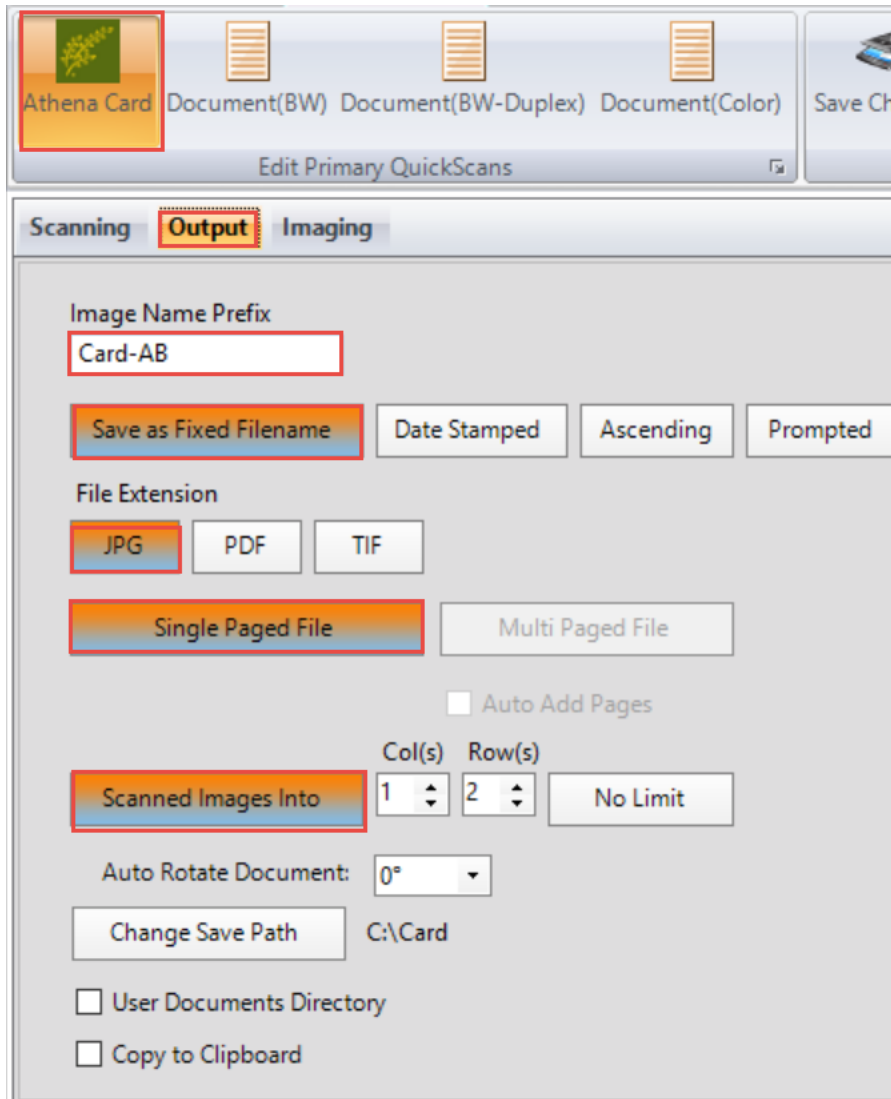


Image Name Prefix: The Card QuickScan default prefix name is Card-AB.

Save as Fixed Filename: All scanned cards will receive the name Card-AB. This means each card scan will overwrite the previous.

File Extension: Card scanning default is JPG.

Single Paged File: Scans images into a single page

Scanned Images Into

Col(s): Number of columns to be displayed in the file.

Row(s): Number of rows to be displayed in the file.


Change Save Path: The default save location for the Card QuickScan is C:\ard.



AmbirScan for athenahealth User Guide

Card – Imaging

The screenshot below displays the default settings for the **Card QuickScan Imaging Menu**:

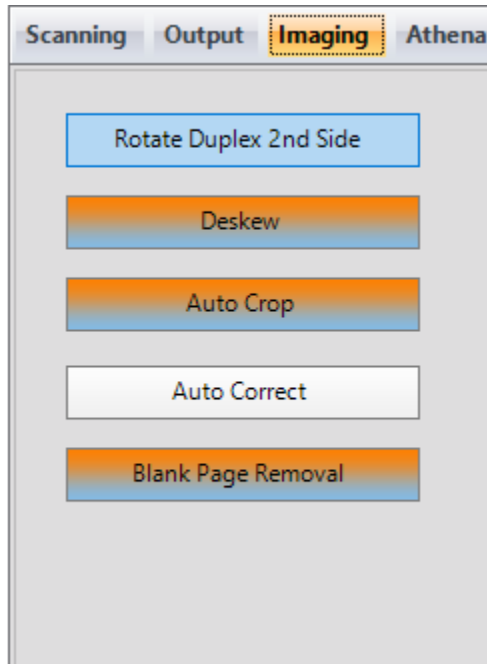

Rotate Duplex 2nd Side:

Deskew: Removes rotation from the scanned media

Auto Crop: Removes black borders from the edge of the scanned media

Auto Correct: Stretches Contrast


Blank Page Removal: Removes blank pages. *This feature is only available for Grayscale and Color scanning.*



AmbirScan for athenahealth User Guide

Document – Scanning

Below is a screenshot of the default settings of the **Document(GS-Duplex) QuickScan Scanning menu:**

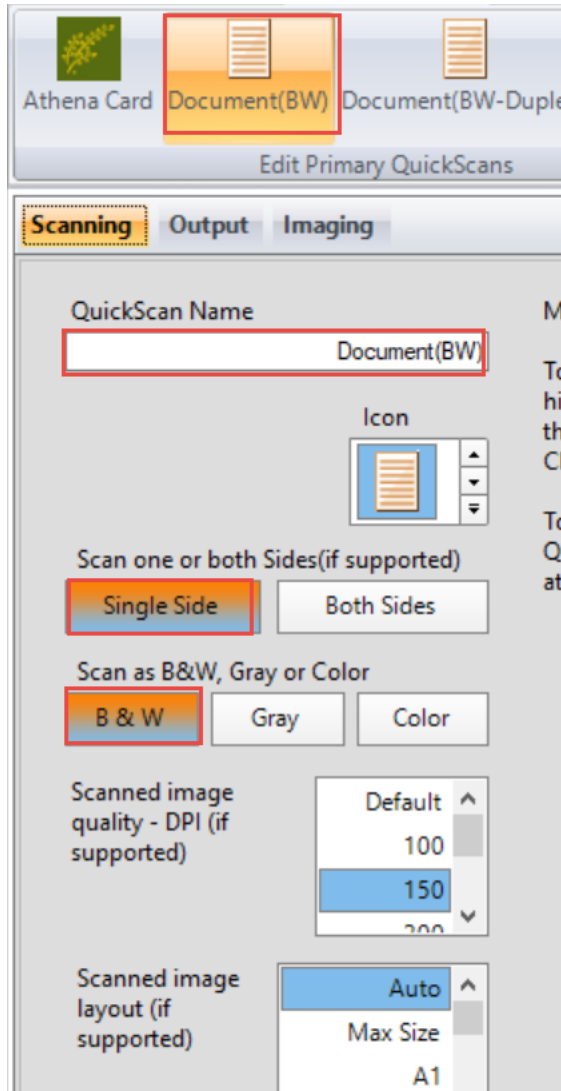

⚠ Please note the settings for each Document QuickScan will vary, depending on their specific functions and color scans.



AmbirScan for athenahealth User Guide

Document – Output

The screenshot below displays the default settings for the **Document(BW-Duplex) Quickscan Output Menu**:

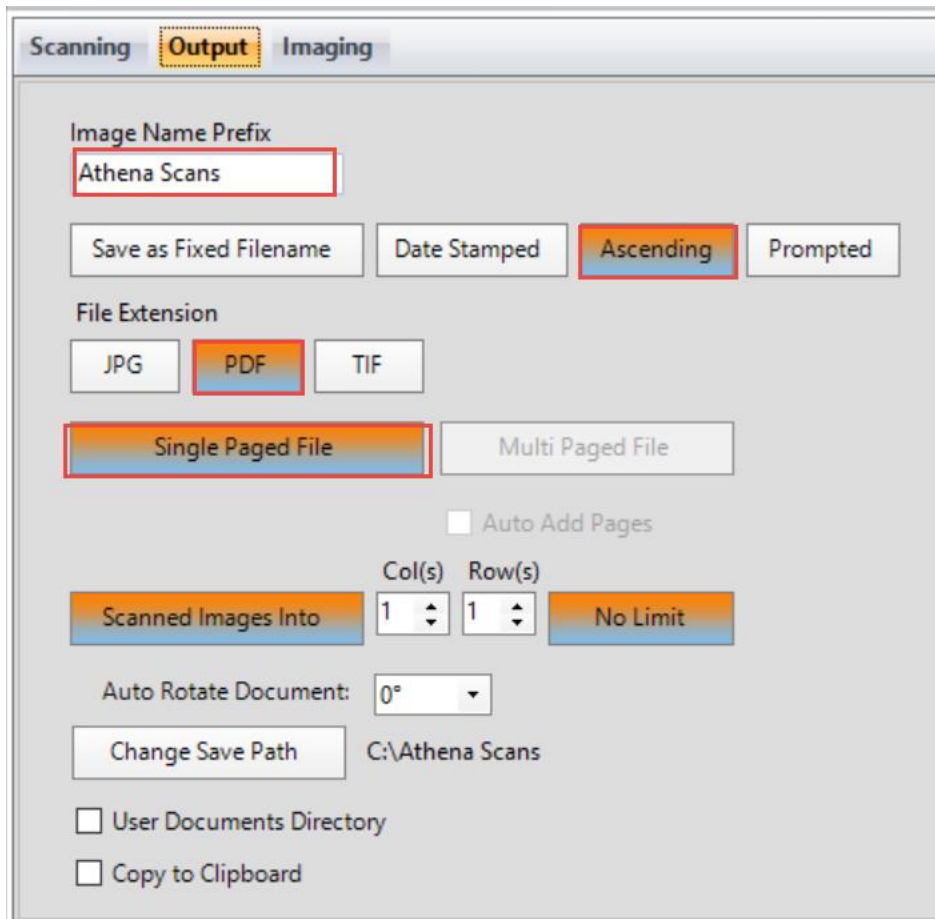


Image Name Prefix: The default prefix for document scans is Athena Scans.

Save as Fixed Filename: Ascending is the default setting for Document QuickScans.

Auto Rotate Document: Adjust the rotation of the document

Single Paged File: Creates a single file for a single page scan

Multi Paged File: Allows scanning of multiple paged documents to be saved in one file

Auto Add Pages: When enabled, allows users to add pages to a multi-page document scan

Change Save Path: The default location for Document QuickScans is C:\Athena Scans. Click on the Change Save Path button to change the default location.

User Documents Directory: Select a path if you wish to save documents locally

Copy to Clipboard: Saves image to a clipboard



AmbirScan®

athenahealth

AmbirScan for athenahealth User Guide

Troubleshooting and Technical Support

Please visit our AmbirScan for Athena page (<http://www.ambir.com/ath-support>) for additional support videos and answers to the most common questions.

For technical support with your AmbirScan for Athena software, please contact Support at (630) 530-5400, option 3 or chat with us online.